



Partial Hand Prosthesis

- ▶ Custom made
- ▶ Lightweight
- ▶ Durable
- ▶ Bendable fingers
- ▶ Four different glove materials

Partial Hand Prosthesis

Features & Benefits

The lightweight Partial hand prosthesis is developed to meet the customers needs of comfort, functionality and natural appearance. It is custom made after the cast of the patients residue limb for maximum comfort and fit. It features bendable fingers to make the hand more functional, the finger stiffness comes in standard, reinforced or super-reinforced.



The partial hand prosthesis is available in four different durable glove materials, PVC, MCV, MCS, and silicone. Each of these materials gives the customer unique features, and combined with the right colour it gives a natural appearance. The gloves can have any colour from the Centri Colour guide, or can be custom coloured from a colour sample of the customer.

See reverse page for order form and ordering information



Centri AB
Kung Hans Väg 2
192 68 Sollentuna
Sweden

T +46 8 505 332 00
F +46 8 505 332 05

www.centri.se
www.fillauer.com
centri@centri.se

CENTRI

Kung Hans Väg 2
192 68 Sollentuna
Sweden

T +46 8 505 332 00 F +46 8 505 332 05
centri@centri.se

ORDER FOR PARTIAL HAND

Customer

Patient

Contact person

Phone

Fax

Finger Stiffness (Thread type)

Standard Reinforced Super-reinforced

Email

Glove type

PVC MCV Silicone MCS Right Left Male Female

Centri Colour Guide or color sample

Veins

Clearly visible (strong) Hardly noticeable (weak)

Length of nails

Medium Long

Zip fastener

No Standard length (from thumb fold) or Length H _____ cm.

A

B

C

D

E

F

G


Note: The most important measurement is from tip of middle finger to distal end of stump (F).

Stiff fingers or other problems?

Remarks:

Continued on back

Make a positive plaster cast of the stump in a **neutral and relaxed position**. Cast should be 3 cm longer than desired glove length. Mark distal end of ulna. We recommend a second cast of the healthy side to show structure.

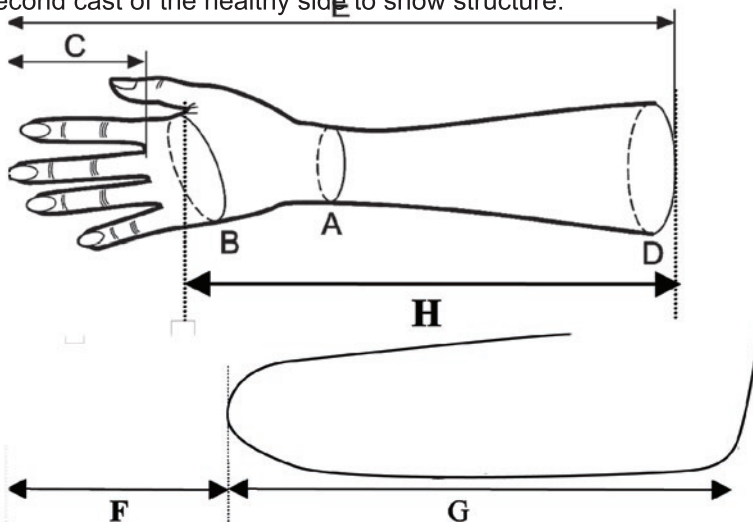

Ordering Numbers

Part no.	Description
613 010 100	Partial Hand with PVC glove - initial order
613 010 200	Partial Hand with PVC glove - repeat order
613 010 500	Partial Hand with MCV glove - initial order
613 010 600	Partial Hand with MCV glove - repeat order
613 010 101	Partial Hand with Silicone glove - initial order
613 010 201	Partial Hand with Silicone glove - repeat order
613 010 501	Partial Hand with MCS glove - initial order
613 010 601	Partial Hand with MCS glove - repeat order